

甘露醇的测定

色谱柱: Hungpu™ Sugar-Ca 300*7.8*8μm


流动相: 去离子纯水

流速: 0.5ml/min

柱温度: 80℃, 检测器温度: 55℃

进样量: 20ul

检测器: 示差折光检测器

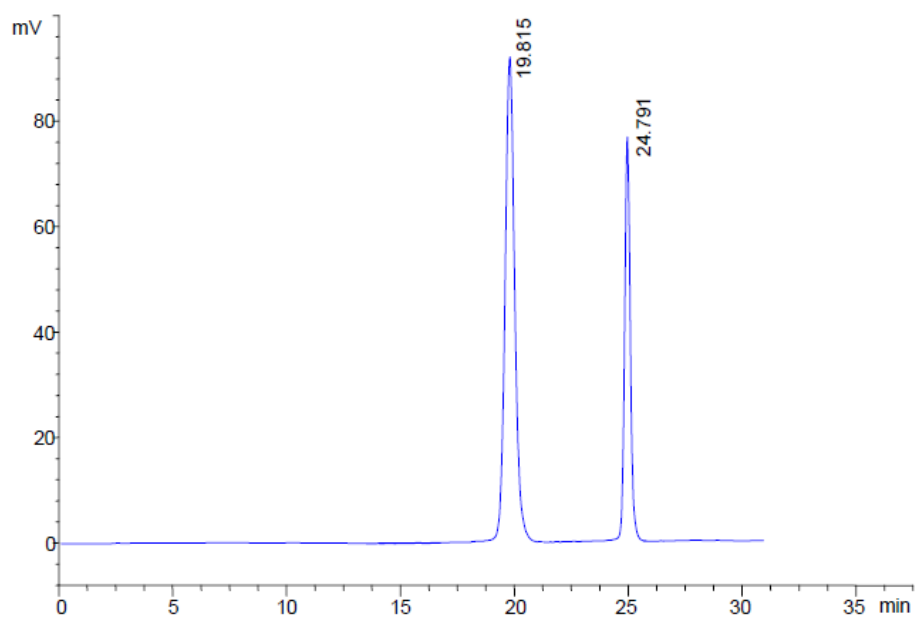

系统适用性实验

Peak #	Ret. Time	Area	Theoretical Plate#	Separation	Tailing Factor
1	19.815	1641577	5193.266	0.000	1.002
2	24.791	1628228	5092.801	3.365	0.994
Total		3269806			

订购信息:

货号	描述
HP-Ca08-7830	Hungpu™ Sugar-Ca 300*7.8*8um

